



Vcheck SAA



For Early Detection of Post-Surgical Infections

The Vcheck Feline SAA test kit is an in vitro diagnostic test kit for the in-clinic, quantitative measurement of Serum Amyloid A (SAA) concentration in feline serum and plasma. In the case of normal felines, the concentration of SAA is very low in their serum and plasma, but it increases when a pathological phenomenon occurs such as infection, tissue damage or inflammatory disorder.


Clinical Applications

- ▶ Monitoring post-operative effects and recovery after surgery
- ▶ Measuring the inflammation response to treatment
- ▶ Detecting the presence of inflammation in cats with no clinical signs

A Closer Look: SAA

The Bionote Vcheck SAA is designed to detect an increased concentration of SAA. This critical detection can indicate and evaluate immunoinflammatory responses in related diseases, infection or tissue damage. This SAA test maintains its high accuracy from previous products, even though the test procedure is much simpler than before.

Diseases that can be monitored with SAA

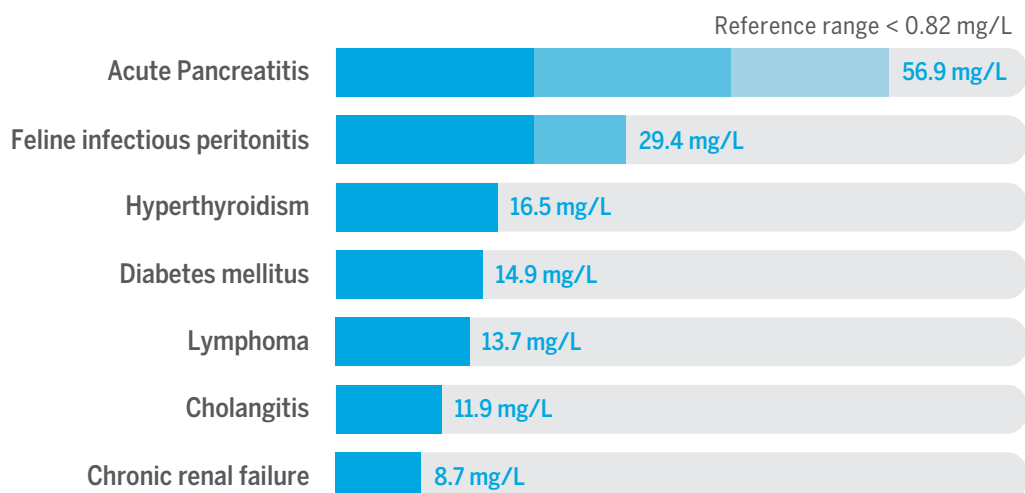

J. Vet. Med. Sci. 70(11): 1247 - 1252, 2008

Specific Clinical Application

The SAA test is critical in many stages of feline care. From regular check-ups to continuous monitoring and post-operative recovery, SAA detection helps diagnose inflammation and infection to provide optimal care for felines.



Regular Check-ups

SAA increases in early stages of inflammation. Test regularly to detect before the presence of symptoms.



Continuous Monitoring

Check SAA after first visit to see if levels have decreased, pointing to a positive response to treatment.



After Surgery

Test before surgery, immediately after surgery and 2-3 times every 12 hours to monitor SAA levels.


Serial measurements are recommended.

*Internal Evaluation Data

Specifications

Species	Feline
Sample Type	Serum, Plasma (heparin) 5 μ l
Measurement	Quantitative
Range	5 - 200 μ g/ml
Testing Time	5 minutes
Storage Condition	1 - 30° C

Simple Testing Procedure

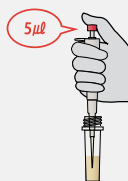

Dilute Sample

Use a 5 μ l pipette to draw 5 μ l of serum or plasma (heparinized) and add into an assay diluent bottle.



Mix

Close the bottle cap and shake 5 - 6 times to mix thoroughly.



Measure

Add 100 μ l of the mixed sample to the sample well of the test device and press [START].

Product Name

Vcheck Feline SAA

Product Number

VCF138DC

Product Type

Device

Packing Unit

10 Tests/Kit



For More Information on
Vcheck V200 or V2400
analyzers visit:
bionote.com
customerservice@bionote.com
800-727-5169

LEARN MORE
about the SAA test

