






Vcheck V200 & V2400

Fluorescent Immunoassay Analyzers for
Quantitative In-clinic Results




Technology

With **Europium**, the stimulating wavelength and emitting wavelength do not overlap, ensuring you receive precise diagnostics for improved care.

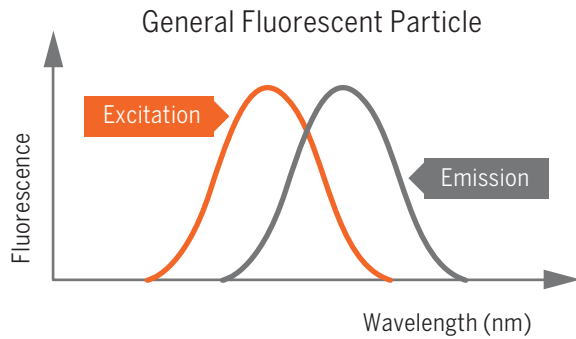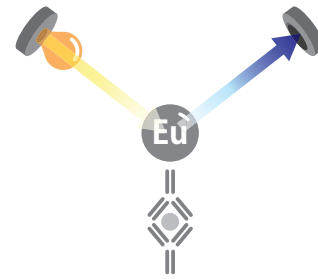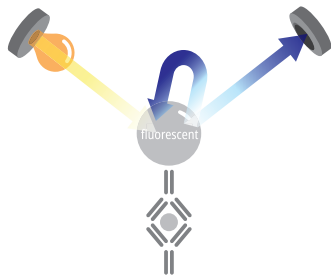

Overlapping Wavelength

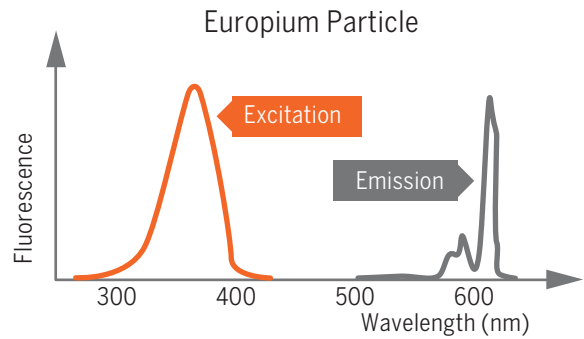

Non-overlapping Wavelength

Vcheck V200



Veterinary Clinics | Hospitals

1 Test at a Time

Dimensions 7.8" wide x 9.4" deep x 8" high

7" Color Touch Screen Display

HLV v2.6(PCD-01) / POCT1-A Connectivity

Built-in Printer



V200 Video

Vcheck V2400



Multi Doctor Practices | Laboratory

24 Tests at Once / 70 Tests per Hour

Dimensions 20" wide x 22" deep x 11.6" high

10" Color Touch Screen Display

HLV v2.6(PCD-01) / POCT1-A Connectivity

Built-in Printer



V2400 Video

Quantitative Tests Available on Vcheck Analyzer

In-House Quantitative Diagnosis Possible

17 Biomarkers Available + More Coming Soon!



Inflammatory
CRP, SAA,
Equine SAA



Pancreatitis
cPL, fPL



Hormone
Cortisol, T4, TSH
Progesterone,
Equine Progesterone



Coagulation
D-dimer



Cardiac
NT-proBNP (Canine, Feline)
TnI (Canine, Feline)



Kidney
SDMA



Antibody
Foal IgG

Clinical Application

Compact and convenient analyzer to expand your in-clinic testing and reduce testing turnaround time.



Regular
Check-Up



Emergency
Check-Up



Pre Anesthetic
Workup



Patient
Prognosis



Combined
Diagnosis

*Pancreatitis: CRP + cPL
*Heart Disease: NT-proBNP + TnI


3-Step Test Procedure



1. Dilute the sample



2. Add diluted sample
into a test device

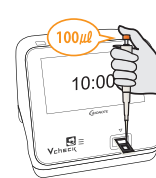

3. Check the test result

Test Reagents and Results



*Soon to be available in North America
Contact your PMI company and request this software interface

Ordering Information

Product Name	Product Number	Product Type	Packing Unit
Vcheck Canine NT-proBNP	VCF132DC	Device	5 Tests/Kit
Vcheck Feline NT-proBNP	VCF130DC	Device	5 Tests/Kit
Vcheck Canine CRP 2.0	VCF109DD	Device	10 Tests/Kit
Vcheck Feline SAA 3.0	VCF138DD	Device	10 Tests/Kit
Vcheck cPL 2.0	VCF129DD	Device	10 Tests/Kit
Vcheck fPL 2.0	VCF127DD	Device	10 Tests/Kit
Vcheck D-dimer	VCF107DD	Device	10 Tests/Kit
Vcheck cCortisol	VCF105DD	Device	10 Tests/Kit
Vcheck T4	VCF106DD	Device	10 Tests/Kit
Vcheck Canine TSH	VCF118DC	Device	5 Tests/Kit
Vcheck cProgesterone	VCF122DD	Device	10 Tests/Kit
Vcheck Canine Tnl	VCF137DC	Device	5 Tests/Kit
Vcheck Feline Tnl	VCF139DC	Device	5 Tests/Kit
Vcheck Equine SAA	VCF141DC	Device	10 Tests/Kit
Vcheck eProgesterone	VCF142DC	Device	5 Tests/Kit
Vcheck Foal IgG	VCF143DC	Device	5 Tests/Kit
Vcheck SDMA	NOW AVAILABLE! VCF125DD	Device	10 Tests/Kit



For More Information on the
Vcheck V200 and V2400
analyzer visit:
bionote.com
customerservice@bionote.com
800-727-5169

LEARN MORE
about the V200 and the
V2400 analyzer

